

R-CU509/R-CU510 Multifunctional Veterinary ICU

RWD



Quiet, Comfortable & Safe



The reference data can be calculated by selecting animal type-age-disease name



Touchscreen LCD operation for best simplicity



Intelligent temperature control system with heating and cooling capabilities



Intelligent Humidity control system with humidification and dehumidification capabilities



Using high-quality oxygen sensor and microcomputer to achieve precise control of O₂



UVC lamp& Anion function for effective sanitation management, to avoid patients cross infection



Automatically clear CO₂ by CO₂ absorbent-internal circulation and ventilation methods, effectively preventing CO₂ accumulation



Built-in failure alarm covering abnormal temperature, humidity and CO₂ levels


From Animal Needs, for Animal Use

Multifunctional Veterinary ICU

Professional Multifunctional Veterinary ICU is indicated for postoperative recovery, shocks, head injury, neurologic disease, tissue hypoxia, cardiopulmonary disease, and numerous other symptoms and situations.

It provides animal patients with the most comfortable medical care environment to help animals recover faster.

Application

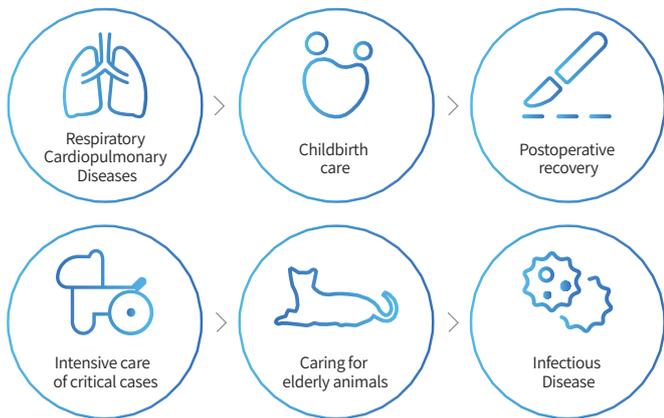

Specifications

Item	Parameters
Temperature Control Adjustment Method	Semiconductor Components
Temperature Control Range	15-40°C(59~104°F), accuracy : ±0.5°C
Humidity Control Range	30-70%RH
O ₂ Concentration Control Range	20-60%
CO ₂ concentration monitoring	500-4000PPM
LED Lamp	Level 1~Level 3
Anion	>8×10 ⁶ PCS/cm ³ (Anion injection port 20cm-30cm)
External Power Outlet	500W
Display	7 inch LCD Touch Screen
Power Supply Voltage	AC 100~120V / 220~240V, 50/60Hz
Weight	50KG (Net weight), 65KG(Packaged weight)
Size	External: W 970 mm * L 530 mm * H 790mm Internal: W 690mm * L 470mm * H 490mm

Order Information

Cat No.	Configuration
R-CU510	Heating/Cooling/ Humidifying/Dehumidifying/ CO ₂ automatic clearing/ oxygen therapy administration
R-CU509	Heating/ Humidifying/Dehumidifying/ CO ₂ automatic clearing/ oxygen therapy administration

RWD Life Science Co.,Ltd

Add: 19-20/F, Building 7A, Shenzhen International Innovation Valley, Dashi 1st Road, Nanshan District, Shenzhen, Guangdong, P.R. China Enquiries: market@rwdstco.com

RWD Life Science Inc.

Add: 850 New Burton Road, Suite 201, Dover, DE 19904 (New Office)
Tel: (858)900-5879 E-mail: service@rwdls.com